Juvora

CERTIFIED  
DISTRIBUTOR

AMANNGIRRBACH


## INSTRUCTIONS FOR USE

CERAMILL PEEK by JUVORA™

CE 1639

2	Dental Prosthesis Forms	33	Norske
3	English	36	Polski
6	Dansk	39	Português
9	Deutsch	42	Româna
12	ΕΛΛΗΝΙΚΑ	45	Slovenčina
15	Español	48	Suomi
18	Français	51	Svenskt
21	Hrvatski	54	Türk
24	Italiano		
27	Magyar		
30	Nederlands		

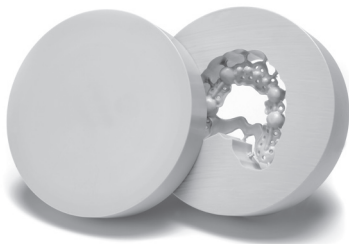1. Telescope

2. Attachment

3. Implant Supported Superstructures

### Indications

---

For the manufacture of frameworks used in the production of:

- i) Full and partial removable dentures and overdentures;
- ii) Implant dentures and implant bars;
- iii) Anterior and posterior crowns, and posterior 3-unit bridges (maximum 1 pontic)

### Contraindications

---

Not to be used in bridges with 2 or more pontics.

### Precautions

---

Follow the milling parameter guidance given in the JUVORA™ processing guide:

- The minimum thickness of secondary crowns and matrix housings should not be less than 0.7 mm
- The minimum thickness for “conventional” retention grids should not be less than 2 mm.

- Minimum size of horizontal retainer in lower jaw 2 mm x 8 mm
- Minimum size of horizontal band in upper jaw 2 mm x 8 mm

### Material Information

---

JUVORA™ dental discs are manufactured from Invibio PEEK polymer

Form:	Solid
Odour:	Odourless.
pH (Value):	Not known.
Boiling Point (°C):	Not known.
Melting Point (°C):	343
Flash Point (°C):	Not known.
Auto Ignition Temperature (°C):	595
Oxidising Properties:	Not applicable.
Vapour Pressure (Pascal):	Not known.
Solubility (Water):	Insoluble.

## Processing Instructions

---

Detailed processing instructions are provided and must be observed by the end users (dental labs and technicians) when milling frameworks from the JUVORA™ dental discs. Dental labs that are supplied with the disc will be certified through Juvora Ltd, a process that involves assessment of their capability to follow the processing guidance to manufacture high quality frameworks and finished dentures where applicable using their existing equipment. This assessment will be completed by Juvora Ltd or their certified milling centres.

Find more information about the processing of PEEK @ [www.juvoradental.com/academy](http://www.juvoradental.com/academy)

## Production

---

JUVORA™ dental discs blanks can only be processed on suitable milling machines with appropriate milling programs. The milling head must be suitable for machining PEEK polymer. Juvora Ltd recommends suitable adaptation of the CAM software by the software manufacturer to allow milling of JUVORA™. For an exact fit, a 5-axis milling machine of stable construction is recommended. Juvora Ltd also recommends cooling the milling head for more accurate milling and to prolong its useful life.\*See foot note.

## Prosthesis removal from the dental disc

---

Remove the prosthesis carefully from the blank using a hard metal cutter, and pay particular attention to the edges of the prosthesis.

### Finishing

---

For finishing work, hard metal cross-cut burs are recommended. Plain cross-cut burs or special fine-toothed burs are particularly suitable. Further fine finishing of horizontal bands and retainers can be carried out with a silicon polisher. Before high-gloss polishing, Juvora Ltd recommends brushing on the polishing motor. The device can be polished to a high gloss with a goat hair or cotton wool buffing wheel. For the insides of telescope attachments, Juvora Ltd recommends fine brushes. Use fine (white) polishing paste.

### Safety Information

---

While machining the JUVORA™ dental discs the following safety precautions are recommended:

- Dust mask or dust extraction
- Personal protective equipment (eye protection, gloves)

### Storage information

---

The Juvora™ dental disc should be stored in dry storage and avoid exposure to direct sunlight. The PEEK polymer from which the Juvora™ dental disc is made is stable and can be stored for an extended period of time (10 year shelf life). It has a working temperature range from cryogenic up to 250°C and hence the storage temperature range for the Juvora™ dental disc is any ambient temperature and humidity.

### Additional Information

---

For additional information contact [info@Juvoradental.com](mailto:info@Juvoradental.com)

\*Note: Follow the instructions of your machine manufacturer regarding the processing of CERAMILL PEEK by JUVORA™. Download processing instructions @ [www.amanngirrbach.com/instruction-manuals](http://www.amanngirrbach.com/instruction-manuals).

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Amann Girschbach AG