

## **MALO CLINIC** study indicates 3X-less bone loss with PEEK-polymer-based dental prosthetics

The MALO CLINIC are taking part in a long term clinical study into the CAD/CAM-millable JUVORA™ disc, made from PEEK-OPTIMA™, as a substructure for the Clinic's signature All-on-4® procedure, where prostheses supported full-arch implant-supported restorations for edentulous patients.



### **THE DENTIST'S PERSPECTIVE**

“Our conclusion, after three years of follow-up, is that it must be the PEEK-OPTIMA™ biomaterial that is helping to enhance patient results, and we think this is really exciting. What JUVORA offers us is confidence that patients are receiving the best possible treatment, and not experiencing unnecessary bone loss. The All-on-4 technique has not changed; our experienced surgeons have not changed – the only change is the material substructure. Therefore, when the JUVORA disc is used as the substructure for All-on-4, the PEEK-OPTIMA material mimics the behaviour of natural periodontal ligament by providing a cushioning and shock-absorbing effect in the mouth, under extreme force.”

Miguel Nobre, Director of Oral Hygiene and R&D, MALO CLINIC

### **THE RESULTS**

At 3 years of patient follow-up, when compared to metal, the JUVORA-based prostheses had positive results<sup>(1)</sup>

- ▶ 3x-lower marginal bone loss
- ▶ Peri-implant disease avoided in most patient cases
- ▶ Patients reported an improved quality of life
- ▶ Absence of correlation between plaque and bleeding
- ▶ High implant survival – close to 100%



## Full-arch implant supported fixed hybrid made with **JUVORA™** with All-on-4® framework data comparison

Data from 3 year clinical study led by MALO CLINIC Lisbon compared with their routine patients receiving a titanium framework implant<sup>(1)</sup>

<b>JUVORA™</b>		<b>Titanium</b>
<b>100%</b>	Implant survival	<b>98%</b>
<b>0.5%</b>	Peri-implant disease (implant level)	<b>2.1-3.0%</b>
<b>0.38mm</b>	Marginal bone loss (maxilla AO4)	<b>1.52mm</b>
<b>1</b>	Bleeding index	<b>2</b>
<b>2</b>	Plaque index	<b>2</b>
<b>9.0</b>	Patient comfort	<b>7.9</b>

### Patient selection

The MALO CLINIC did not discriminate in its patient selection and included a number of cases with complex histories, including heavy bruxers. This means that even in complex cases, JUVORA exhibits better performance than conventional dental materials such as titanium.

# Handle extreme forces with JUVORA™

## Shock-absorbing properties of JUVORA

PEEK-OPTIMA™ offers several mechanical benefits over metal frameworks, namely, a more favorable strength-to-weight ratio and the potential for shock absorption. Metal prosthetics are naturally stronger; however they are less resilient and can be less comfortable for the patient<sup>[2]</sup>. If left untreated, bone loss can lead to pain, infection and loss of lower facial structure as the infection eats away at the jaw bone.

The shock-absorbing properties of PEEK-OPTIMA are explained by its modulus of elasticity, which is very similar to that of cortical bone. When stress transfer is measured against conventional dental materials, PEEK exhibits the most shock absorption<sup>[3]</sup>.



[1] Nobre, M; Guedes, M; Silva A. European Association for Osseointegration (EAO) Congress Symposium: What if dental prosthetic could influence quality of life? September 2019. [2] Kurtz, S (2012). "PEEK Biomaterials Handbook". Edition 1, Oxford, Elsevier Inc. [3] Based on the tendency of a material to deform when opposing forces are applied as expressed by the average elastic (Young's) modulus values of the two materials (Titanium 110 GPa or PEEK 4GPa). \* All-on-4 is a registered trademark of Nobel Biocare.

# Hear it from the MALO CLINIC


Click buttons below to learn more:



Antonio Silva, Master Technician, MALO CLINIC discusses how the material properties of PEEK-OPTIMA™ enable JUVORA™ to create a new future of dental prosthetics



Dr Miguel Nobre, Director of Oral Hygiene and R&D, MALO CLINIC, presents a lecture on the clinical evidence of PEEK in dentistry

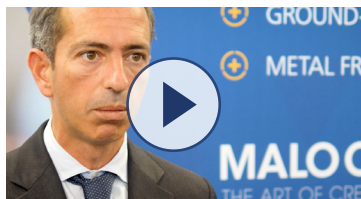

Dr Carlos Moura Guedes, MALO CLINIC, Clinical Director explores why materials hold the key to the evolution of the Clinic's signature All-on-4™ treatment protocol

Victrex plc and/or its group companies ("Victrex plc") believes that the information contained in this document is an accurate description of the typical characteristics and/or uses of the product or products, but it is the customer's responsibility to thoroughly test the product in each specific application to determine its performance, efficacy, and safety for each end-use product, device or other application. Suggestions of uses should not be taken as inducements to infringe any particular patent. The information and data contained herein are based on information we believe reliable. Mention of a product in this document is not a guarantee of availability.

Victrex plc reserves the right to modify products, specifications and/or packaging as part of a continuous program of product development. Victrex plc makes no warranties, express or implied, including, without limitation, a warranty of fitness for a particular purpose or of intellectual property non-infringement, including, but not limited to patent non-infringement, which are expressly disclaimed, whether express or implied, in fact or by law.

Further, Victrex plc makes no warranty to your customers or agents, and has not authorized anyone to make any representation or warranty other than as provided above. Victrex plc shall in no event be liable for any general, indirect, special, consequential, punitive, incidental or similar damages, including without limitation, damages for harm to business, lost profits or lost savings, even if Victrex has been advised of the possibility of such damages regardless of the form of action.

Supporting information is available on request for all claims referenced in this document.

Copyright © 2020 Invisio Ltd. INVIBIO™, PEEK-OPTIMA™, INVIBIO BIOMATERIAL SOLUTIONS™ are trademarks of Victrex plc or its group companies. All rights reserved.